BEAUTIFLY

IllumiEyes Therapy

Blush

series

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

[/beautifly.eu](https://www.facebook.com/beautifly.eu) [@beautifly.eu](https://www.instagram.com/beautifly.eu) [beautifly.eu](https://www.beautifly.eu)

[/beautiflyeu](https://www.youtube.com/channel/UC...) [@beautifly.eu](https://www.tiktok.com/@beautifly.eu)

f1

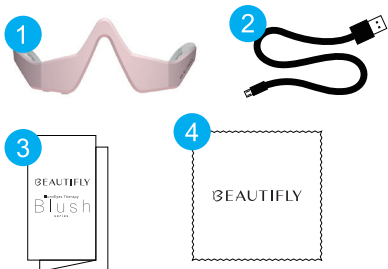f2

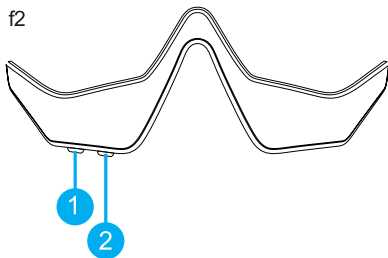f3

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Size:	139.6*53.5*37mm
Weight:	42g
Frequency:	0.5V
Charging voltage:	5V
Charging current:	120mA
Maximum operating current:	120mA
EMS frequency:	60-1 kHz

Set includes: [see page 2, f1]

1. EMS eye massage device
2. USB cable
3. User manual
4. Cleaning cloth

Device parts description: [see page 2, f2]

1. On/Off button /switching operating mode
2. Intensity change button

Before use

Check the battery by long pressing the on/off button, if the device does not turn on, connect the USB charging cable.

How to use?

ON-OFF button/mode switch

Power on: Long press the ON/OFF button for 2 seconds.

Power off: long press the ON/OFF button for 2 seconds.

Mode switching: Click the power button to switch to the first mode. Click the button again to go to the second mode. Another press activates the first (default) mode.

1. Turn on the device and put on **like glasses**, the first working mode is EMS long red light - continuous wave. Press the button again to switch to pulse wave mode.
2. There is a time limit of 10 seconds in the working state. If no operation is performed within 10 seconds, the device will automatically turn off.
3. Charging: The device is charging by connecting with a USB cable, the red light is on when charging. When fully charged, the light shines on green.

The operation of the continuous wave mode lasts 10 minutes after this time, the device will turn off.

Pulse wave operation lasts 10 minutes after this time, the device will turn off.

Precautions

- It is normal that there is a slight tingling sensation when using the device, the intensity depends on the moisture of the skin, the higher the humidity, the stronger the feeling.
- Do not use this device for purposes other than those specified in the user manual.
- Do not use the device for people with acute, chronic diseases, fever, heart disease, pregnancy, facial surgery, skin rash, blood diseases, menstruation, people with pacemakers, orthodontic appliances.
- Do not use the device around damaged skin, mouth or private parts of the body.
- Protect the device from impacts.
- Keep the device away from dust, vibration, hot or too cold places.
- If the packaging is damaged, make sure that the device is working properly before use.
- Keep the device out of the reach of children.
- In case of failure, do not try to disassemble the device yourself.
- Do not submerge the device in water, do not hold it under a stream of water - the device is not waterproof.
- Use only the USB cable supplied with the product.

Storage

- Store the device in a dry and clean place.
- Keep away from high temperatures.
- Keep out of reach of children.

Cleaning

- Wipe with a damp cloth.
- Do not submerge the device in water - it is not waterproof.

More information about the product is available on the manufacturer's website at: [see page 2, f3]

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EN: Electrical and electronic devices, as well as batteries and accumulators, must not be mixed with other household waste. This is indicated by the crossed-out trash can symbol on the packaging. Used electrical and electronic devices, batteries, and accumulators must be taken to local selective waste collection points or to the retailer, and more information can be found in your municipality. Electrical and electronic equipment, batteries, and accumulators may contain hazardous substances, mixtures, and components harmful to the environment and human health. Improper handling or damage to them can lead to health hazards or environmental

pollution during further disposal or recycling.

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EN: The manufacturer hereby declares, at its sole responsibility, that the LumiEyes Therapy Blush product complies with the Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, as well as with the Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU with regard to the list of restricted substances. The full text of the declaration of conformity is available at the following website:

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